



Molift Transfer Pro

Enhancing seated transfers - empowering mobility

Molift Transfer Pro

A transfer is more than just a movement - it is an important part of the patient's rehabilitation and maintenance of physical abilities. Molift Transfer Pro contributes to a safe and healthy care environment for both the patient and the caregiver.

Optimized for seated transfers

Molift Transfer Pro is a sit-to-stand device for seated transfers, e.g. between bed, wheelchair, toilet or similar. The user should be able to follow instructions, pull up their weight and maintain good core stability.

Moreover, it is an ideal solution for various healthcare settings such as hospitals, nursing homes and home care.

Enhanced comfort and safety by innovative design

Featuring an ergonomic and spacious handle that allows multiple grip options for both patient and caregiver, and easy access to the basin for washing hands, Molift Transfer Pro ensures enhanced usability. The stable seat pads provide a comfortable and safe transfer while the soft padded leg support adapts to the shins, increasing user comfort.

The user can easily place their feet on the high-friction footplate thanks to the low step-on height and sloped edge. For caregivers, adjusting the width of the legs is easy, and the design enables the device to come close to the seating area, facilitating a safe transfer. The central brake offers both auditory and visual feedback, securing both wheels simultaneously.

On the side of the handle are attachment points for the Molift Assist accessories.



Handle

Multiple ergonomic grip options for the patient and caregiver, ensuring both comfort and ergonomic positioning of hands.



Central brake

Intuitive one-foot touch central brake, locking both wheels, regardless of their positioning, while also providing both auditory and visual feedback.



Seat pads

The stable seat pads enhance a comfortable transfer for the user. The incut edges provide a good grip, making folding them up and down easy for the caregiver.



Turning radius

High quality wheels for exceptionally easy manoeuvring meaning minimal force is needed when pushing/pulling the device.

Clinical benefits:

- Reduces work-related injuries or strain on the caregiver.
- Enables early rehabilitation.
- Comfortable and safe seated transfers.
- Supports the user to use their physical ability.
- Enables single-handed care.

Natural movement pattern

A natural movement pattern refers to the way humans naturally move their bodies during activities in daily life.

These natural movement patterns are inherent to human physiology and have been shaped by evolution to facilitate activities such as standing up and sitting down.

Molift Transfer Pro supports the natural movement pattern by allowing feet to be placed under the user, and the spacious design of the handle makes it possible for the user to lean forward when getting up from the seat.



Scan QR code to learn more

Product	Item no.	Specifications										
Molift Transfer Pro	E16090300	<table border="0"> <tr> <td>Height</td> <td>46.5 inches</td> </tr> <tr> <td>Length</td> <td>36.25 inches</td> </tr> <tr> <td>Width</td> <td>26.25 inches</td> </tr> <tr> <td>Foot plate height</td> <td>2.25 inches</td> </tr> <tr> <td>Turning radius with open/closed legs</td> <td>40.5/45 inches</td> </tr> </table>	Height	46.5 inches	Length	36.25 inches	Width	26.25 inches	Foot plate height	2.25 inches	Turning radius with open/closed legs	40.5/45 inches
Height	46.5 inches											
Length	36.25 inches											
Width	26.25 inches											
Foot plate height	2.25 inches											
Turning radius with open/closed legs	40.5/45 inches											
User height Minimum - Maximum	4 feet 11 inches to 6 feet 6.75 inches											
Weight Total weight Maximum user weight (SWL)	62.25lbs. 440lbs.											
Material Powder coated aluminum, plastic, steel, polyuretan foam		Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution.										




Accessories



Molift Soft Handles 2 pcs

Helps the user reach the handle.

Molift Soft handles.....E82511

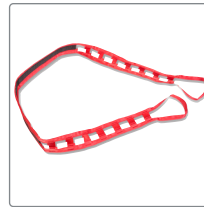

Molift Assist Strap

A safety strap for extra support and safety. Available in two sizes, S/M and L/XL.


Molift Assist Strap S/M.....N82508
Molift Assist Strap L/XL.....N82509



Molift Assist Sleeve

An Assist Strap with sliding sleeve, to assist user to rise up/sit down and support them during transfer.

Molift Assist Sleeve S/M.....N82528
Molift Assist Sleeve L/XL.....N82529

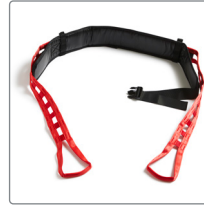

Molift Assist Belt

An assist Strap with a sliding sleeve with buckle to assist the user to rise up/sit down and support them during transfer.

Molift Assist Belt S/M.....E82964
Molift Assist Belt L/XL.....E82965



Molift Assist Sleeve Wide

An Assist Strap L/XL with firm sliding sleeve to give more support to assist user to rise up/sit down and support them during transfer.

Molift Assist Sleeve WideE90301



Molift Transfer Pro IV-Pole

Height-adjustable IV-Pole to be mounted on the Molift Transfer Pro.

Molift Transfer Pro IV-PoleE1630030



Molift Transfer Pro O-Bracket

The oxygen tank bracket is for securely storing and transporting oxygen tanks.

Molift Transfer Pro O-BracketE1630040

Etac is a world-leading developer of innovative assistive devices and patient handling equipment. Since 1973, we have been committed to improving quality of life for the individual, family members and caregivers.

For the intended purpose of the products and latest information, please visit www.hovermatt.com/molift-transfer-pro